

**PHILIPS**  
RESPIRONICS

Sami the Seal

Reassuringly  
**reliable**

Sami the Seal combines a child friendly compressor design, with our highly efficient SideStream nebulizer, to provide fast treatments to children and help support aerosol therapy compliance.


**Robust and reliable**

The system has been designed to cope with the demands of delivering frequent treatments in the hospital, and is also a robust and reliable option for the home.

**Key features**

- Includes Tucker the Turtle character mask
- Partnered with proven SideStream technology
- Robust design suitable for use in hospital and home



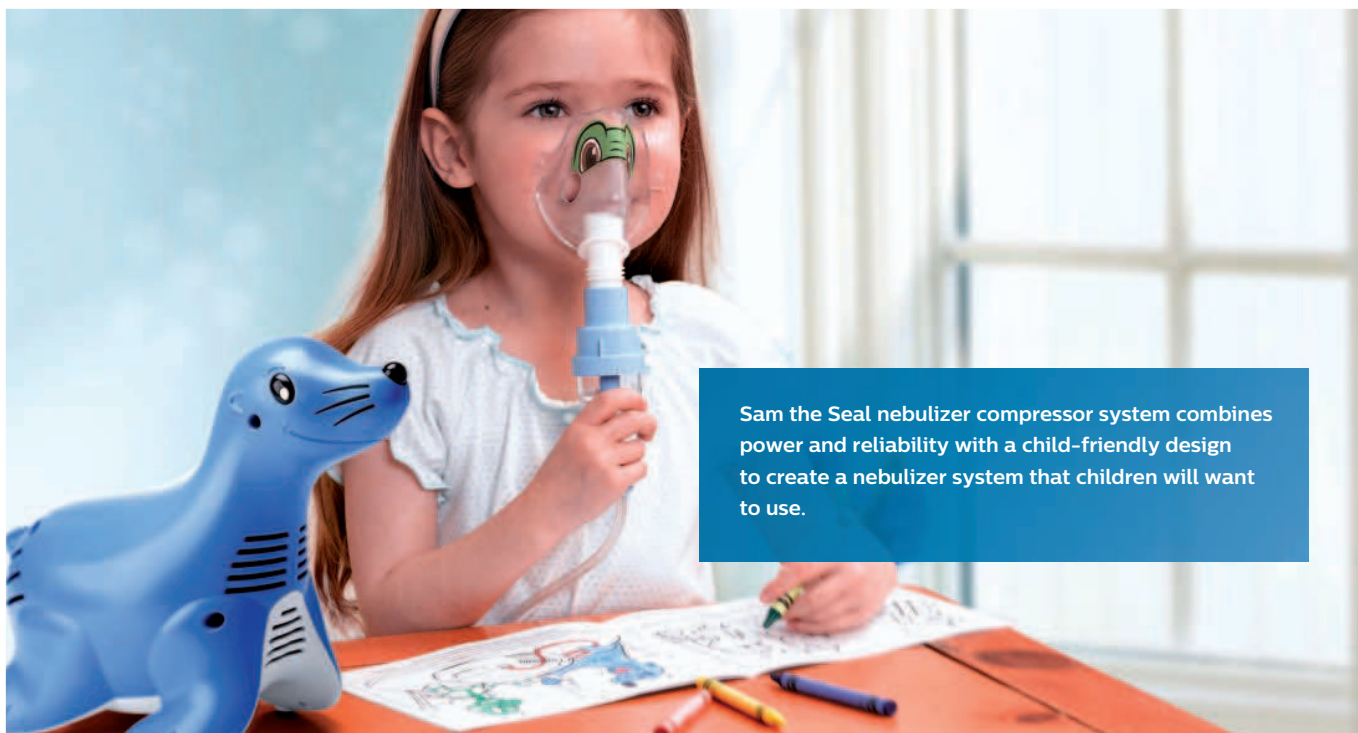

Sam the Seal nebulizer compressor system combines power and reliability with a child-friendly design to create a nebulizer system that children will want to use.

## Cute without compromise

Sami is rated for continuous use, making it suitable for frequent use in hospital or in the home. Sami also includes a built-in thermal protector safety feature.

### SideStream reusable nebulizer

SideStream can be used with commonly prescribed respiratory medications and offers fast and efficient drug delivery. SideStream has an active venturi that provides air flow in addition to that supplied by the compressor, resulting in short treatment times.\*

### Tucker the Turtle pediatric mask

Our pediatric character mask has been designed to help reduce aerosol deposition in and around the eyes.\*



Product description (with voltage/frequency)	Part
Sami the Seal compressor 230 V 50 Hz Euro plug	1093268
Sami the Seal compressor 220 V 60 Hz UK plug	1093270
Sami the Seal compressor 220 V 60 Hz Korean plug	1093237
Sami the Seal compressor 220 V 60 Hz US plug	1093235
Reusable pediatric replacement nebulizer kit	1101484

### Product specification

Thermal overload protector	Resettable thermal fuse; functioning temperature 130°C
Performance flow	6 lpm @ 10 psi
Maximum pressure	29~43.5 psi
Size (l x w x h)	267 mm x 155 mm x 226 mm
Weight	1.6 kg
Operating mode	Continuous
Warranty	3 years

© 2015 Koninklijke Philips N.V. All rights are reserved. Philips Healthcare reserves the right to make changes in specifications and/or to discontinue any product at any time without notice or obligation and will not be liable for any consequences resulting from the use of this publication.



Philips Respironics  
1010 Murry Ridge Lane  
Murrysville, PA 15668 USA

Philips Respironics Asia Pacific  
+65 6882 5282

Philips Respironics Australia  
+61 (2) 9947 0440  
1300 766 488 (toll free, Australia only)

Philips Respironics United Kingdom  
+44 800 1300 845